
SDI

POLA LIGHT

ADVANCED TOOTH
WHITENING SYSTEM

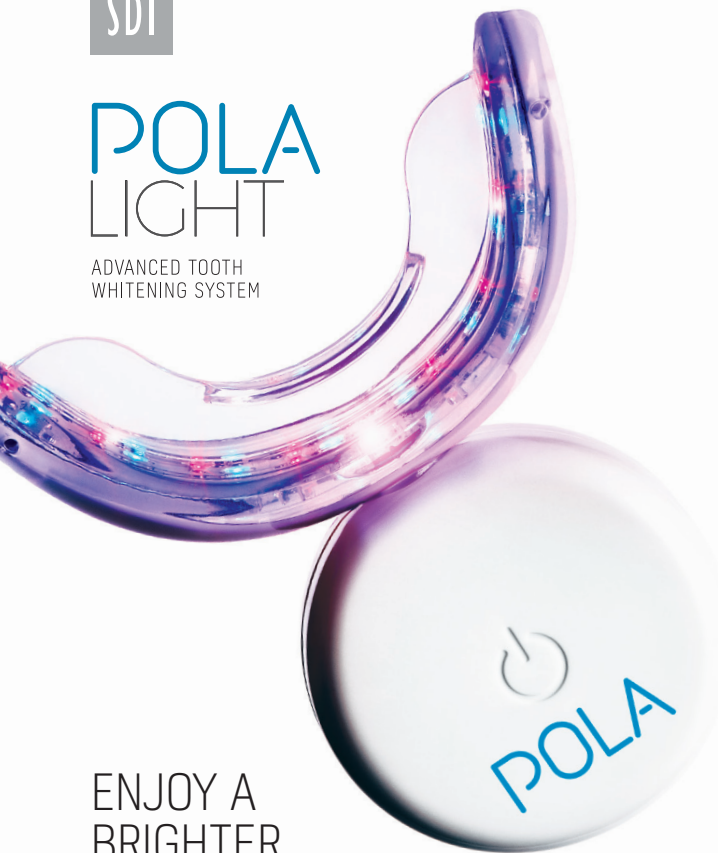ENJOY A
BRIGHTER,
MORE CONFIDENT
SMILE IN 5 DAYS

EASY
TO USE

FROM JUST
20 MINUTES
A DAY

WHITER
TEETH IN
FIVE DAYS

PROFESSIONALLY
USED AND
RECOMMENDED

PROFESSIONAL TOOTH WHITENING SYSTEM

Millions of people across 100 countries have trusted Pola to achieve a safer, quicker and whiter smile.

The Pola Light advanced tooth whitening system combines Pola's award winning whitening formula, with targeted LED technology. And all from the comfort of your home.

- Fast acting, delivering you whiter teeth in 5 days
- Formulated to safely remove long term stains
- Returns your smile to your whitest natural shade possible
- Extra hydration to ensure maximised patient comfort
- Remineralization to strengthen your teeth
- Fully rechargeable LED light for your future tooth whitening treatments

Inclusions

1 x Pola LED mouthpiece

4 x 3g [0.1 oz.] Pola Day or Pola Night tooth whitening syringes

1 x USB charging cord


1 x Instruction guide

DIRECTIONS FOR USE

- 1** Clean your teeth and rinse your Pola LED mouthpiece before use.

Optional: Record your tooth shade on the last page.

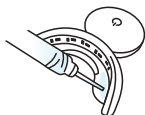- 2** Attach the syringe tip provided to the syringe. Apply a small line of Pola Day or Pola Night gel evenly along the inner front of the top and bottom of your Pola LED mouthpiece.

Using the markings on the syringe as a guide, allow one fifth of the syringe (0.6g) per tooth whitening session on both top and bottom of the tray.

- 3** Fit the Pola LED mouthpiece to your mouth and bite down.

Refer to the label on your Pola Day or Pola Night whitening syringe for your whitening formula type, then use the following duration instructions:


Pola Day 9.5% hydrogen peroxide whitening

- Press 'on' button twice for the **Blue/Red** LED
- LED will activate for 10 minutes

Pola Night 22% carbamide peroxide whitening

- Press 'on' button once for the **Blue** LED
- LED will activate for 15 minutes

When the LED activation ends, your whitening session is complete. Repeat twice a day for 5 days.

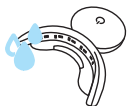- 4** After reaching your treatment duration, remove the mouthpiece and wipe away any excess gel/saliva with a dry tissue. Rinse the mouthpiece and wash it with lukewarm water. Do not submerge it in water.

- 5** Rinse mouth or brush teeth once whitening session is complete.

FREQUENTLY ASKED QUESTIONS

HOW MUCH WHITER WILL MY TEETH BE?

Many teeth whitening manufacturers make claims of '5 shades' whiter or similar. However, this can be misleading as it depends on the original shade of your teeth and the current staining on your teeth. Pola is formulated to remove the toughest stains from your teeth and return your teeth to their whitest natural color.

HOW LONG WILL MY TEETH STAY WHITE?

Your lifestyle is the greatest influence on your oral health. If you are an avid drinker of coffee, wine and/or soda, your teeth whitening will diminish more quickly than someone who prefers water. Typically, your whitening will be at its best for the first 6 months, before slowly diminishing.

HOW MUCH POLA DAY OR POLA NIGHT WHITENING GEL DO I APPLY?

Only a small amount is required for each session. Do not overfill your mouthpiece, as gel may spill on to your gums and cause sensitivity. Using the markings on the syringe as a guide, allow one fifth of the syringe (0.6g / 0.02 oz) per tooth whitening session on both top and bottom of the tray.

THE TREATMENT GUIDELINES RECOMMEND TWO APPLICATIONS PER DAY. SHOULD I COMPLETE BOTH APPLICATIONS IN ONE SITTING?

For best results, we recommend completing the two applications one after the other.

IS POLA WHITENING SAFE?

Yes. Clinical studies have shown that teeth whitening with carbamide peroxide or hydrogen peroxide is safe if performed under the guidance and supervision of a trained dentist. With a history of improving teeth for 50 years, our Pola whitening scientists have refined an industry leading tooth whitening formula that removes tough 'rusted on' stains, while remaining completely safe for your teeth.

WILL THE WHITENING PROCESS HURT?

A small minority of people may experience short term teeth sensitivity during the whitening process. This sensitivity will subside quickly and does not lead to tooth or gum problems. If you are experiencing pain, stop the whitening process and consult your dentist.

HOW DO I CHARGE THE POLA LED MOUTHPIECE?

Once the LED light starts to flash, charge the Pola mouthpiece using the supplied USB charger cord. The light will continue to flash while charging and will stop flashing once fully charged.

HOW SHOULD I STORE THE POLA LIGHT KIT?

For best results, store the Pola whitening gel in a safe location with a temperature between 35°F - 77°F [2°C - 25°C].

HOW CAN I ORDER MORE POLA WHITENING GEL?

To order more Pola whitening gel for future treatments, simply contact your dentist. Your dentist can provide you with additional Pola whitening gel syringes, which are compatible with the Pola Light system.

HOW DO I CLEAN THE POLA LED MOUTHPIECE?

Please use only water and soap to clean the Pola mouthpiece. Avoid the use of abrasive materials and cleaning substances.

INSTRUCTIONS FOR USE FOR PATIENTS AND DENTAL PROFESSIONALS

Pola Day 9.5% hydrogen peroxide Pola Night 22% carbamide peroxide Advanced tooth whitening system

Pola Day and Pola Night are neutral pH take home tooth whitening solutions. Each is pleasantly flavoured, crystal clear with a unique blend of soothers, conditioners and a water content high enough to ensure maximised patient comfort.

INDICATIONS FOR USE:

- Whitening of discolored vital teeth.
- Whitening of discolored non-vital teeth.

HYDROGEN PEROXIDE EQUIVALENT

- Pola Night 22% Carbamide Peroxide = 7.3% Hydrogen Peroxide

STORAGE AND SHELF LIFE

- Store in a cool place and out of direct sunlight at temperatures between [2°-25°C / 35°-77°F].
- Do not use after expiration date.
- A warranty period of one (1) year is provided for the Pola LED mouthpiece. Please retain proof of purchase for you to claim the warranty.

FOR EFFECTIVE RESULTS

- Do not eat, drink or smoke during treatment.
- Do not smoke immediately after treatment - wait for at least two hours.
- Foods and drinks containing strong colors should be avoided for at least 48 hours after the treatment period or consumed in moderation.
- Use gel at room temperature.

PRECAUTIONS

- To be used under the supervision of a dentist.
- Keep out of reach of children and pets.
- Not to be used by pregnant or lactating women.

- Patients with a history of chemical allergies are advised to carry out allergy testing by a specialist before using these products.
- Discontinue use if any unusual sensitivity or reactions occur and consult your dentist.
- Do not inject Pola Day or Pola Night into the body.
- CAUTION: FEDERAL LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A DENTIST.
- SDS available at www.sdi.com.au or contact your regional representative.

Note: Patients with sensitivity will require a day or two break between treatments. Shortening treatment time is also suggested.

WARNING FOR CANADA ONLY

- If irritation (such as redness, swelling, soreness) of the gums or the mouth occurs, discontinue use and consult an oral health professional.
- Products containing peroxides are not recommended for use by children under 12 years of age.
- Consult your oral health professional before prolonged use of this product.
- Avoid swallowing the cosmetic or part thereof.
- Avoid contact of the product with the eye.
- Avoid direct contact of the active surface of the tooth whitening product with the gums and/or salivary flow.

FIRST AID

Eye (contact): Hold eyelids apart and flush with running water for at least 15 minutes. Seek medical advice if symptoms persist.

Skin (contact): Wash thoroughly with water and seek medical advice if symptoms persist.

Ingestion (large amount): If swallowed, do NOT induce vomiting. Rinse mouth. Drink plenty of milk or water and seek medical advice.

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SDI Limited