

## REQUIRED MATERIALS:


- Teeth model – trim base of model so that it slightly extends past the gingival border
- SDI Tray material
- SDI Spacer Resin
- Curing light
- Vacuum former

## PROCEDURE

1. Apply approximately 0.5 mm of Spacer Resin to the teeth on the model selected for treatment. Apply to about 2.5 - 3.5 mm from the gingival line.



2. Light cure each tooth for 20 seconds.

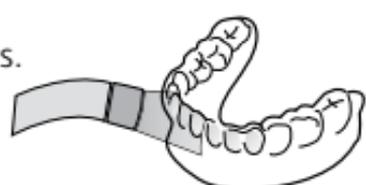


3. Firmly attach a piece of tray material into the holder. Note: We recommend the tray holder to be fitted with a rubber gasket. This will prevent slipping of the material during heating.

4. Position the tray holder for heating.

5. Place the model on the centre of the base.

6. Position the pre-heated heating element over the tray.

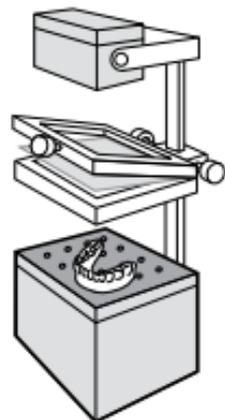


7. Heat until the material starts to sag. Turn the heat off. Note: The thickness of the tray will depend on a number of factors, e.g. heating time, position of heating element above tray etc.

8. Lower the tray onto the model.


9. Apply vacuum for approximately 10 seconds.

10. Allow the model to cool before removing from vacuum former.

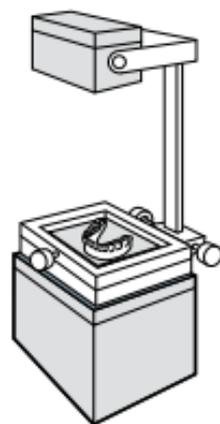

11. Carefully remove the tray away from the impression.

12. Trim the tray to 0.5mm from the gingival margin.

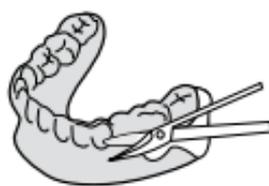

**CAUTION: Ensure that the dispensing tip hub of the Spacer Resin is firmly attached to the syringe by twisting the hub securely onto the syringe.**

Precautions: Avoid prolonged contact of Spacer Resin with skin when dispensing, as it may cause inflammation or skin sensitization. Any persons having known resin allergies should immediately discontinue the use of Space Resin. MSDS for Spacer Resin is available at [www.sdi.com.au](http://www.sdi.com.au) or contact your regional representative.

NOTE: On first usage, or after prolonged storage, extrude a small amount of Spacer Resin onto a mixing pad for familiarity with the viscosity and rate of extrusion. Use at room temperature. Do not expose material to direct light.

**Spacer Resin 2 syringe refill**

Reorder 7700012  
2 x 1g Spacer Resin syringes  
accessories

**Tray material - 1mm refill**

Reorder 7750005  
20 piece refill