SDI polaprofessional 35%

ADVANCED ENDODONTIC WHITENING SYSTEM

INSTRUCTIONS FOR USE

E
N
G
L
I
S
H

Pola Professional 35% is a hydrogen peroxide whitening system designed to medically bleach endodontically treated teeth.

INDICATIONS:

- Whitening of discoloured non-vital teeth that have been endodontically treated.

COMPOSITION: Contains 35% Hydrogen peroxide.

NET WEIGHT:

Depending on kit purchased. Please refer to outside of box:

- 0.3g Bleach powder
- 2.25g Hydrogen peroxide liquid
- 1g Gingival Barrier

SUGGESTIONS FOR SAFER NON-VITAL BLEACHING:

1. Verify adequate endodontic obturation.

The root canal filling must have been completed to appropriate endodontic standards. The root canal should be completed at least a month before bleaching. The quality of root canal obturation should always be assessed clinically and radiographically prior to bleaching. Adequate obturation ensures a better overall prognosis of the treated tooth. It provides an additional barrier against damage by peroxide to the periodontal ligament and periapical tissues.

2. Use protective barrier.

This is essential to prevent leakage of bleaching agents which infiltrate between the gutta-percha and root canal walls, reaching the periodontal ligament via dentinal tubules, lateral canals or the root apex. The ideal barrier should protect the dentinal tubules. **Caution:** Seal gutta-percha to a level which is coronal to potential communication with periodontal tissue. We recommend using a glass ionomer cement such as Riva SC (SDI Ltd).

3. Avoid acid etching.

Acid etching has not proven to be beneficial (Casey et al 1989) and the use of caustic chemicals in the pulp chamber is undesirable, as periodontal ligament irritation may result.

4. Recall patients periodically.

Bleached teeth should be frequently examined both clinically and radiographically. Root resorption may occasionally be detected as early as 6 months after bleaching. Early detection improves the prognosis as corrective therapy may still be applied.

TIPS:

- Isolate tooth effectively. Use of a rubber dam is advised.
- Protect oral mucosa to prevent chemical burns.
- Advise the patient that a whitened non-vital tooth will not look the same as a vital tooth.
- Delay any esthetic restorations for 2 to 4 weeks after treatment.
- Advise patient not to smoke, consume highly coloured, hot or cold foods or beverages for at least 2 days.
- Recall patients periodically to examine the treated teeth (clinically & radiographically).
- Large diameter suction tips are not recommended as they can spread the gel to other areas.
- A patient consent form is available for download at www.sdi.com.au

PROCEDURE:

- Take a pre-operative photograph to record tooth shade.
- Take pre-operative X-ray to verify the presence of an acceptable root canal treatment and the absence of apical pathology.
- Isolate tooth/teeth undergoing whitening using a rubber dam. Ensure the dam fits tightly at the cervical margin.
- Open lingual access to pulp chamber and remove all restorative material and ensure correct endodontic access has been achieved. Remove gutta-percha to 1-1.5 mm below the cemento-enamel junction (CEJ).
- Place a 1 mm thick protective barrier eg a conventional glass ionomer over the exposed gutta-percha.
- Take a periapical radiograph to check that the barrier has been well placed.
- Clean pulp chamber with water.
- Open powder pot. Take one Pola Professional 35% syringe, **firmly** attach a tip, and **carefully pull back** plunger to release pressure.
- Carefully** extrude contents of syringe into the pot.
- Immediately mix using a brush applicator until gel is homogeneous.
- Place the Pola Professional 35% gel into the opened pulp chamber.
- Leave gel in for 8 min (Optional: curing light may be used at this stage, please refer to light manufacturer's recommendations) NOTE: If using a curing light, two curing lights can be used (one on each side of tooth) to activate Pola Professional 35%.
- Suction gel out.
- Repeat steps 11 to 13 for up to 4 applications. Note: One pot should only be used for 2 applications or within 20 minutes of mixing.
- Coat the inner labial surface of access cavity with bleach gel, place dry cotton pellet. Seal with composite material (we recommend Glacier (SDI Ltd) in a contrasting shade) without using a bonding agent (this helps prevent excess pressure build up within tooth).
- Verify occlusal clearance.
- Recall patient 2-3 days later for assessment. Procedure can be repeated until desired shade is achieved - first place rubber dam, remove temporary cover and gel then repeat steps 15 and 16 using freshly mixed gel. Note: The longer the tooth has been discoloured, the longer it will take for the bleaching treatment to remove the discoloration.
- Remove gel from pulp chamber and wash with water to remove residue.
- Place a temporary seal. (we recommend a glass ionomer cement such as Riva Protect (SDI Ltd) in a pink shade).
- Wait 10-14 days before completing final restoration.
- Record final shade.

STORAGE AND HANDLING:

- Store in cool place (2°-25°C / 35°-77°F).
- Do not use after expiration date.

PRECAUTIONS:

- For professional in-office use only.
- Patient must wear protective eyewear.
- Operators must wear gloves, mask and protective eyewear.
- Do not use on pregnant or lactating women.
- Do not use on children under 14 years of age.
- Do not use the Gingival Barrier on any persons having known resin allergies.
- Do not use Pola Professional 35% on any persons having known peroxides allergies.
- Anyone with a history of chemical allergies is advised to undergo allergy checks before treatment.
- Pola Professional 35% Liquid syringe may be under pressure – remove cap carefully.
- Do not allow Pola Professional 35% mixed gel or liquid to come in contact with skin, eyes and soft tissues.
- Pola Professional 35% will not lighten any restorative materials.

FIRST AID:

- SKIN/TISSUES (contact):** Immediately wash thoroughly with water. Apply a slurry of sodium bicarbonate and water to the affected area. Seek medical attention if symptoms persist.
- EYES (contact):** Open the eye wide and thoroughly wash for 15 minutes with running water. Seek medical attention if symptoms persist.
- INGESTION:** Rinse mouth with water. Gargle with salt water and drink lots of milk or water. Seek medical attention if symptoms persist.
- INHALATION (powder):** Remove to fresh air. If breathing is difficult, seek medical attention.

PRODUCT RANGE:

Pola Professional 35% 1 Patient Kit
2 x 0.3g Pola Professional 35% powder pots
2 x 2mL Pola Professional 35% syringes
1 x 1g Gingival Barrier syringe
Accessories


Pola Professional 35% 3 Patient Kit
6 x 0.3g Pola Professional 35% powder pots
6 x 2mL Pola Professional 35% syringes
3 x 1g Gingival Barrier syringes
Accessories

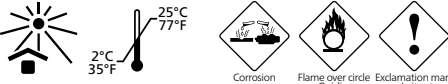Made in Australia by SDI Limited
Bayswater, Victoria 3153
Australia 1 800 337 003
Austria 00800 0225 5734
Brazil 0800 770 1735
France 00800 0225 5734
Germany 0800 1005759
Ireland 00800 0225 5734
Italy 00800 0225 5734
New Zealand 0800 734 034
Spain 00800 0225 5734
United Kingdom 00800 0225 5734
USA & Canada 1 800 228 5166
www.sdi.com.au
www.polawhite.com.au

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