



AURA EASYFLOW

UNIVERSAL FLOWABLE COMPOSITE



4

EASY SHADES

CONTROLLED
FLOW,
HIGH
STRENGTH

CONTROLLED FLOW, HIGH STRENGTH

Aura Easyflow is a light cured flowable composite, refined to provide the right features for every layer of use. Whether it is used as a radiopaque liner under direct restorations or superficially to repair defects in esthetic zones, Aura Easyflow is engineered to maximize clinical success.

ENGINEERED FILLER TECHNOLOGY

Aura Easyflow is based on the successful filler technology utilised in SDI's Aura range of composites. Its state of the art filler system defines its versatility for multiple clinical needs to ensure high strength and natural esthetics.

The particle size distribution has been tailored by SDI scientists to polish to a high gloss that is maintained over time. The hybrid combination improves mechanical properties without affecting the required flow¹.

COMPRESSIVE
STRENGTH

421MPa


4
EASY SHADES



CONTROLLED FLOW, DOES NOT SLUMP

Producing conservative restorations that preserve the tooth, often leads to structural irregularities on cavity walls. A flowable composite must adapt into every irregularity to form a complete seal at the tooth-restorative interface. Filler loading affects the flowability of a composite¹ and at 56% [filler load by weight], Aura Easyflow is a superb choice to achieve precise flow around the cavity walls.

Aura Easyflow's tip design for syringes, makes it easy to inject the material into cavities of any depth and any location within the mouth. Tactically, the clinician can be confident that the material is placed with precision.

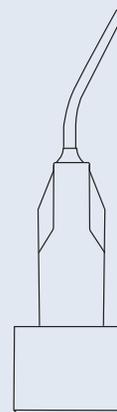

HIGH RADIOCAPACITY

HIGH VISIBILITY ON RADIOGRAPHS

Over time, it is important to evaluate the contours of restoration and distinguish between the restorative material and secondary caries¹. Aura Easyflow contains barium glass, selected for its high radiopacity and therefore, high visibility on radiographs.

Aura Easyflow makes it easier to see the distinction between tooth and restorative material, enabling confident clinical diagnosis throughout the life of the restoration.

**AURA
EASYFLOW**

265%_{Al}


ENAMEL³

220%_{Al}

DENTIN³

120%_{Al}

RADIOCAPACITY*

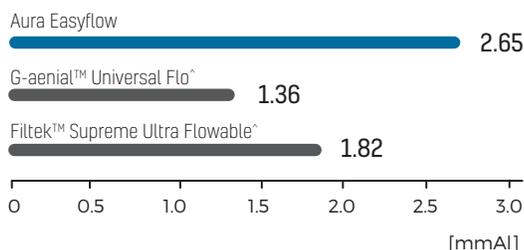


» In Class V restorations, composites with lower modulus are desired as they are capable of flexing during tooth function which may reduce stresses along the bonding agent interface and the likelihood of debonding. «

STRONG MECHANICAL PROPERTIES


Aura Easyflow's strength asserts its use as a versatile restorative material. With an impressive compressive strength of 421MPa, Aura Easyflow's ability to withstand forces of mastication, meets that of conventional composites and is superior in the flowable composite category.

Flexural strength influences the use, stability and clinical success of a restoration². Aura Easyflow displays sufficient flexural strength. The material has enough flexural strength to be used on the surface of Class I and Class II restorations that are often under masticatory stress, while maintaining enough flexibility to be used in Class V scenarios.


COMPRESSIVE STRENGTH*

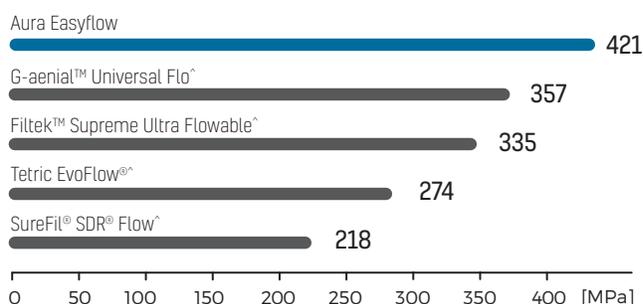

FLEXURAL STRENGTH*

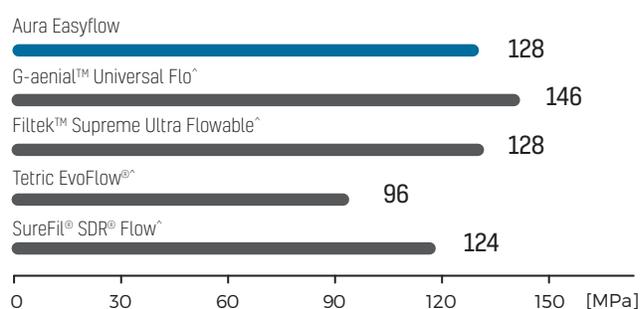

FLEXURAL MODULUS*

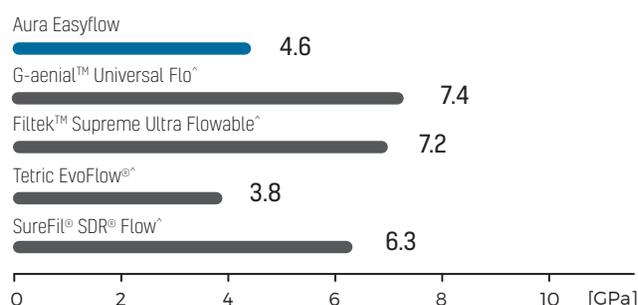


EASY SHADE MATCHING

EASY SHADE MATCHING

Aura Easyflow is designed with four unique shades to simplify colour selection in esthetic zones and enable better inventory management. Shade selection is simplified by three methods:

1. Shades are arranged according to chroma, from weaker intensity to higher intensity of color.
2. Shades are equidistantly spaced, creating visual logic to the eye.
3. Shades have a single opacity, making it easier to predict the final esthetic of the restoration.

Four unique Aura Easyflow shades cover 16 VITA® A1-D4® shades, enabling an easy switch from traditional shading systems.

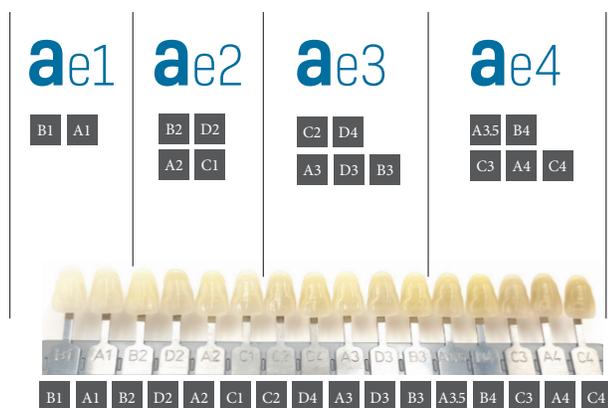

Approximate equivalence to VITA® A1-D4® shade guide

LIFELIKE OPTICAL PROPERTIES

Aura Easyflow contains the same optical properties as SDI's Aura, Aura Easy and Aura Bulk Fill. The development of Aura Easyflow completes the Aura composite range to provide a restorative material for every direct clinical use.

Aura Easyflow's filler and resin technology is specially tailored to display a chameleon effect. There is a lifelike blend of translucency, opacity, opalescence and fluorescence to mimic that of an unrestored, natural tooth.

KEY FEATURES

4 shades cover 16 VITA® A1-D4® shades

Higher radiopacity than enamel

Controlled flow, does not slump

Flexural strength compatible for Class V

Very high compressive strength

Engineered filler technology

Easy shade matching

Lifelike optical properties

Stable esthetics over time

INSTRUCTIONS

1 Clean and isolate tooth



2 Etch tooth surface with Super Etch 37% phosphoric acid for 20 seconds



3 Remove excess water. Keep moist



4 Apply bonding agent to saturate all internal surfaces according to manufacturer's instructions



5 Directly inject the flowable composite in increments of 2mm or less in: Class V restorations, pit & fissure sealants, conservative Class I, II, III and IV restorations or other indications as required



7 Cure Aura Easyflow for 20 seconds in increments of 2mm




8 Polish and finish



STABLE ESTHETICS OVER TIME

Flowable composites are often selected as a material of choice for the repair of superficial defects. Particularly in the esthetic zone, it is important that the shade remains stable for the life of the restoration. Aura Easyflow demonstrates color stability over time, which facilitates long lasting satisfaction with the esthetics of the restorative work undertaken.


AURA EASYFLOW COLOR STABILITY: COLOR CHANGE (ΔE) OVER TIME FOR A DISC STORED IN DI WATER AT 37°C*

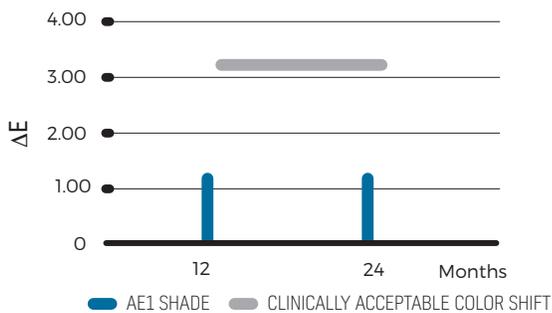

TECHNICAL TABLE

Filler load (total)	56% weight 32% volume
Depth of cure (mm)	2
Compressive strength (MPa at 24 hours)	421
Flexural strength (MPa at 24 hours)	128
Volumetric shrinkage (%)	4.8
Radiopacity (%Al)	265

ORDER DETAILS

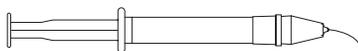

SYRINGES

Aura Easyflow Syringe Refill

1 x 2g syringe
5 x single use disposable tips

Ae1	8566010
Ae2	8566011
Ae3	8566012
Ae4	8566013

TIPS

Flowable Disposable Tips 20 x single use disposable tips (20 gauge)	8100160
---	---------

Flowable Disposable Tips 200 x single use disposable tips (20 gauge)	8100161
--	---------

* SDI internal data

^ G-aenial™ Universal Flo, Filtek™ Supreme Ultra Flowable, Tetric EvoFlow® and SureFit® SDR® Flow are not the registered trademarks of SDI

1 Kusai Baroudi et al (2015), 'Flowable Resin Composites: A Systematic Review and Clinical Considerations', Journal of Diagnostic Research, NCBI

2 Sayna Shamszadeh et al (2013), 'Comparison of Flexural Strength of Several Composite Resins Available in Iran', Journal of Dental School

3 KM Lachowski et al (2013), 'Study of the radiopacity of base and other liner dental materials using digital radiography system, DentoMaxillioFacial Radiology, NCBI

~ Yap AU et al (2002), 'Changes in flexural properties of composite restoratives after gaining in water', Journal of Operative Dentistry, 27:468-74



YOUR SMILE. OUR VISION.



02-2024

M100428 C



+D036M1004281S

MADE IN AUSTRALIA
by SDI Limited
Bayswater, Victoria 3153
Australia
www.sdi.com.au

AUSTRALIA 1800 337 003
AUSTRIA 00800 0225 5734
BRAZIL 0800 770 1735
FRANCE 00800 0225 5734
GERMANY 0800 100 5759

ITALY 00800 0225 5734
NEW ZEALAND 0800 808 855
SPAIN 00800 0225 5734
UNITED KINGDOM 00800 0225 5734
USA & CANADA 1 800 228 5166