

SDI | ULTRAMAT 2

HIGH SPEED CAPSULE MIXER



HIGH SPEED, MULTI-USE CAPSULE MIXER

The Ultramat 2 is a capsule mixer. With high speed, it has a mixing frequency of 4600 oscillations per minute. The 5 year warranty assures your purchase is a long lasting investment.

MULTI-USE MIXER

The multi-use Ultramat 2 is suitable for most brands of composite, glass ionomer, amalgam and luting cement capsules.

SAFETY FEATURE

The micro-switch automatically halts trituration when the lid of the Ultramat 2 is lifted.

SIMPLE TOUCH CONTROL SYSTEM

QUIET

EASY TO CLEAN



FEATURES

VARIABLE TIMER

The Ultramat 2 features a 1-16 second microprocessor-crystal controlled oscillator timer.

MIXING FREQUENCY

The Ultramat 2 has a mixing frequency of 4600 oscillations per minute +/-2%.

DUAL VOLTAGE

The Ultramat 2 operates on voltages between 120V~ and 240V~, 50Hz and 60Hz, for compatibility with multiple power voltages and frequencies.

1-16
SECOND
MICROPROCESSOR

4600
OSCILLATIONS

**120 &
240V~**

ADVANTAGES

5 year limited warranty

1-16 second adjustable timer

Quiet

For mixing most leading encapsulated dental materials

Safe and simple to use

TECHNICAL SPECIFICATIONS

Width	217mm	8.5"
Depth	173mm	6.75"
Height	180mm	7.12"
Voltage	120/240V~	50/60Hz
Weight	2.765kg	6.95lbs

5 YEARS
WARRANTY




CAPSULE COMPATIBILITY

The Ultramat 2 can be used for the trituration of encapsulated dental restorative materials from SDI Limited as well as from other companies. The following chart lists some examples of SDI encapsulated dental restorative materials.

DENTAL DEVICE (CAPSULES)	TRITURATION TIMES (SECONDS)
1, 2, 3, 5 spill Lojic+	6 ± 1
1, 2, 3 spill Permite / 1, 2, 3 spill GS-80	8 ± 1
1, 2, 3, 5 spill Ultracaps+	8
5 spill Permite / 5 spill GS-80	6 ± 1
Ultracaps S	9
Riva Light Cure / Riva Light Cure HV / Riva Bond LC / Riva Protect / Set / Riva Luting / Riva Luting Plus / Riva Self Cure / Riva Self Cure HV / Stela Capsule	10



ORDER DETAILS

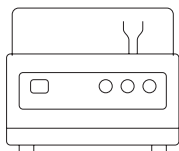

EQUIPMENT

Ultramat 2 (USA - includes power cord for USA)	5546045
Ultramat 2 (All countries - not including power cord for the region the products is shipped to, user to source power cord with equivalent electrical rating or above)	5546052



YOUR SMILE. OUR VISION.



ACN 008 075 581



MADE IN AUSTRALIA
by SDI Limited
Bayswater, Victoria 3153
Australia
www.sdi.com.au

AUSTRALIA 1800 337 003
AUSTRIA 00800 0225 5734
BRAZIL 0800 770 1735
FRANCE 00800 0225 5734
GERMANY 0800 100 5759

ITALY 00800 0225 5734
SPAIN 00800 0225 5734
UNITED KINGDOM 00800 0225 5734
USA & CANADA 1 800 228 5166