

# ultrammat S

HIGH SPEED AMALGAMATOR

AMALGAMADOR DE ALTA VELOCIDADE  
AMALGAMADOR DE ALTA VELOCIDAD  
WYSOKOOBROTOWY MIESZALNIK  
HOCHFREQUENZ MISCHGERÄT  
VIBREUR ELECTRONIQUE A GRANDE VITESSE  
VIBRATORE PER AMALGAMA AD ALTA VELOCITÀ  
HOGE TOEREN CAPSULESCHUDDER  
AVANCERAD HÖGHASTIGHETS KAPSELBLANDARE  
HØJHASTIGHEDS AMALGAMATOR  
NOPEA VAIKUTTEINEN AMALGOIJA  
AMALGAMRYSTER MED HØYENERGI  
ΥΨΗΛΗΣ ΤΑΧΥΤΗΤΑΣ ΔΟΝΗΤΗΣ ΑΜΑΤΑΜΑΤΟΣ  
YÜKSEK HIZLI AMALGAMATÖR  
VYSOKORYCHLOSTNÍ TŘEPAČKA KAPSLÍ  
HITRI MEŠALEC AMALGAMA  
ВИСОКО СКОРОСТНА АМАЛГАМОБЪРКАЧКА  
VYSOKO RÝCHLOSTNÝ AMALGAMÁTOR  
NAGY SEBESSÉGÜ KAPSZULARÁZÓGÉP  
高速アマルガムミキサ-  
高速银汞调拌机



SDI

S stands for simple!

# ultramats



## High speed amalgamator


- Accurate and consistent trituration
- Low vibration and extremely quiet
- Simple to use - just 3 buttons
- Suitable for all your capsules
- Easy to clean
- Low power indicator and shut off function
- Dual voltage



# ultrammat S features:


## So simple, just 3 simple buttons

There are three time setting options, 6, 8 and 10 seconds to mix all your encapsulated materials. The unit can be stopped at any time by pressing these buttons during operation or just by lifting the lid.

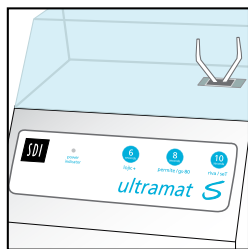

## So easy to clean!

The hygienic one surface cover facilitates infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.

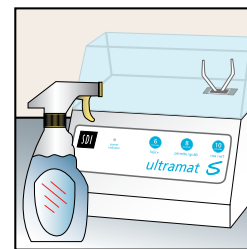


## Accurate and consistent trituration

The Ultramat S features a microprocessor controlled oscillator timer for accurate and consistent trituration.




## Dual voltage

The Ultramat S operates on voltages between 100-120V and 220-240V, 50 and 60Hz for compatibility with multiple power voltages and frequencies.

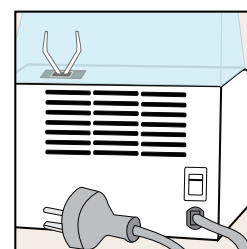


## Low vibration and extremely quiet

Unlike some amalgamators, the Ultramat S has a low noise output. Minimal vibration of all external parts ensures it will not move from its standing position.

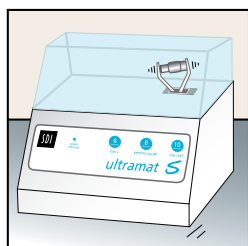


## Low power indicator and shut off function

Should the power supply fall outside the required ranges of 100-120V and 220-240V, the power indicator will flash blue and the unit will not function until correct power is supplied to it. This prevents under-mixed materials.

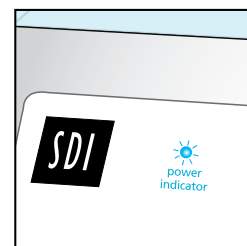


## Suitable for all your capsules

The bendable forks can fit a multitude of round end capsules (amalgam capsules), and also those with protruding nozzles (glass ionomer cement capsules).

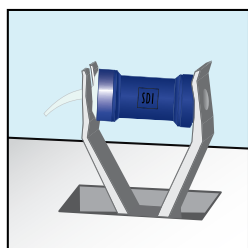



## Safety feature

The Ultramat S will immediately stop mixing when the transparent lid is opened.

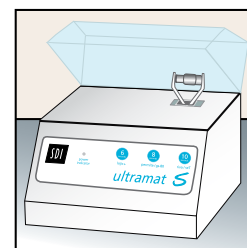

## High energy output and low energy input

Possessing a mixing frequency of 4550 oscillations per minute (+/- 1%), the Ultramat S ensures complete mixing of any dental material, however it consumes relatively low power.

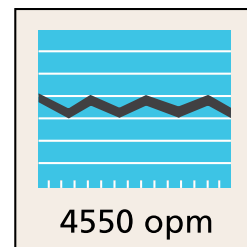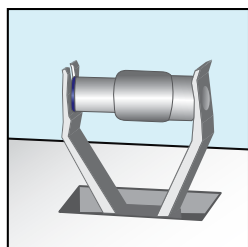

# ultramats



#### Dimensions:

width	depth	height
217mm	173mm	180mm
8 1/2"	6 3/4"	7 1/8"

#### Voltage:

100 - 120 volts

220 - 240 volts

50 - 60 Hz

#### Weight:

2.43kg

5.4lbs


#### Ultramat S

(All countries - includes power cord  
for region product is shipped to)

Reorder 5546058



Made in Australia by SDI Limited  
Bayswater, Victoria 3153  
Australia 1 800 337 003  
Austria 00800 2255 734  
Brazil 0800 770 1735  
France 00800 022 55 734  
Germany 0800 100 5759  
Ireland 01 886 9577  
Italy 800 780625  
Spain 00800 022 55 734  
United Kingdom 00800 022 55 734  
USA & Canada 1 800 228 5166  
[www.sdi.com.au](http://www.sdi.com.au)

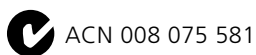

+D036M55462013

M554620 C