



# AURA BULK FILL

UNIVERSAL BULK FILL COMPOSITE



FILL,  
PACK  
AND GO




# THE SIMPLICITY OF BULK FILLING

## SPEEDY AND SIMPLE POSTERIOR COMPOSITE RESTORATIONS

- Natural-looking finish with just one universal shade.
- High depth of cure without the need for a capping layer.
- Easy placement in a matter of seconds.
- The waxy consistency is packable and non-sticky to optimise handling.


### FILL, PACK AND GO

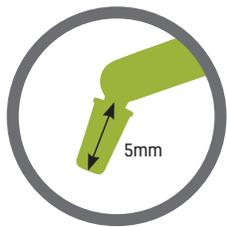

ONE BULK FILL  
PLACEMENT UP  
TO 5mm

20sec  
LED LIGHT CURE

5mm  
BULK FILL

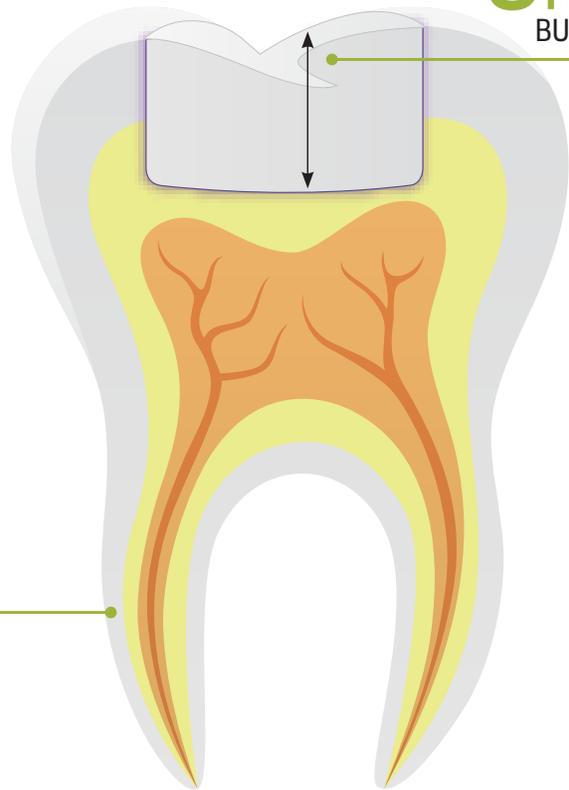

## ACHIEVABLE AESTHETICS CASE STUDIES

DR BOB LOWE (USA)



DR. GREGOR THOMAS (GERMANY)



# ONE SHADE. ONE FILL.

## LOW VOLUMETRIC SHRINKAGE\*

Aura Bulk Fill delivers reduced stress to provide long-lasting restorations.

## EXCELLENT FLEXURAL STRENGTH\*

Aura Bulk Fill resists major deformations without fracturing.

## HIGH COMPRESSIVE STRENGTH\*

Aura Bulk Fill resists the high forces of compression in the mouth.

2.21%

115MPa

375MPa

\*SDI internal data

## INDICATIONS

Posterior restorations

Core buildup

Splinting


Sandwich technique with glass ionomer





## DEPTH OF CURE

The opacity of Aura Bulk Fill is a function of the difference in the refractive index of its two components, the filler and the resin.

The curing process alters the refractive index of the resin marginally, to match the refractive index of the filler. This lowers the opacity temporarily, allowing deeper light penetration for a high depth of cure. Upon curing, the indices move apart again to arrive at the ideal opacity for Aura Bulk Fill's brilliant chameleon-like finish.

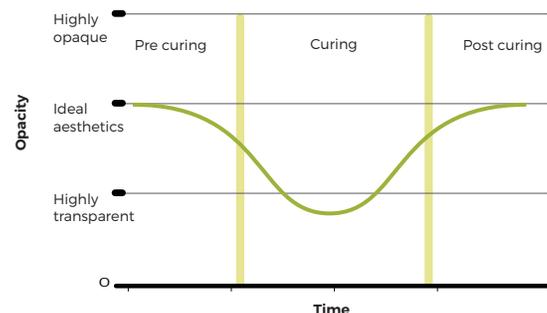

# ORDER DETAILS

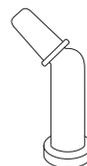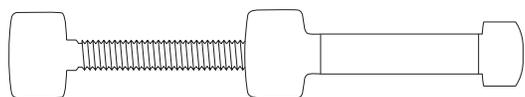

## SYRINGES

Syringe Refill - 4g

8560012

## COMPLETS

Compleat Refill - 20x 0.25g = 5g total

8565012

» I have been using Aura Bulk Fill for three years now. The resin comes out smooth with no pull back and cures extremely strong. The polish is amazing and my posterior composites have never been better. The characteristics are ideal for any restoration, primary or permanent teeth. I have recommended this product to many associates and I plan to continue using the Aura Bulk Fill as my go to composite. «

**DR LIEL GRINBAUM** DMD

## AURA PRODUCT RANGE

Aura  
Aura Easy  
Aura Bulk Fill



YOUR SMILE. OUR VISION.



MADE IN AUSTRALIA  
by SDI Limited  
Bayswater, Victoria 3153  
Australia  
[www.sdi.com.au](http://www.sdi.com.au)

**AUSTRALIA** 1800 337 003  
**AUSTRIA** 00800 0225 5734  
**BRAZIL** 0800 770 1735  
**FRANCE** 00800 0225 5734  
**GERMANY** 0800 100 5759

**ITALY** 00800 0225 5734  
**NEW ZEALAND** 0800 808 855  
**SPAIN** 00800 0225 5734  
**UNITED KINGDOM** 00800 0225 5734  
**USA & CANADA** 1 800 228 5166