

COMPLET APPLICATOR

PREMIUM / ECONOMICAL

CENTRIX
C-R® Syringe System
Precision Placement System for all Dental Materials

INSTRUCTIONS FOR USE

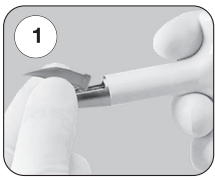

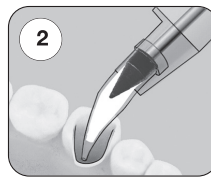

ENGLISH

CENTRIX C-R SYRINGE SYSTEM,
SDI COMPLET APPLICATOR PREMIUM/ECONOMICAL
For precision placement of dental materials.

1) Place into the syringe barrel.
2) Squeeze syringe with a slow, steady pressure.
3) Remove and dispose after use.

Do not use chemi-clave on plastic parts. Do not soak or wipe disposable and plugs in alcohol or other solvents.

CLEANING

TYPE	PROCEDURE
Ultrasonic	Remove all parts and ultrasonic for 15-20 minutes at 122°F (50°C)
Manual	Hand-scrub using a soft bristle in a mild detergent. Rinse under sterilize water and stand vertically to air dry before sterilizing.

STERILIZATION

TYPE	PROCEDURE
Cold solution	Prepare solution following manufacturer's instruction. Immerse the applicator in the solution for length of time as recommended by the cold solution manufacturer. Dry as per instructed.
Steam	Wrapped: 121°C for 15 min OR 132 °C for 10 min



Manufactured by
CENTRIX INCORPORATED
770 River Road Shelton, CT 06484 USA
(203) 929-5582, (800) 235-5862
FAX (203) 929-6804
www.centrixdental.com

EC REP
CENTRIX EUROPE
Hauptstrasse 64 B,
65719 Hofheim,
Germany

EMERGO AUSTRALIA
Level 20, Tower II, Darling Park, 201 Sussex Street
Sydney, NSW 2000
Australia

Distributed by
SDI Limited
3-15 Brunsdon Street
Bayswater, Victoria 3153
Australia

Australia 1 800 337 003
Austria 00800 0225 5734
Brazil 0800 770 1735
France 00800 0225 5734
Germany 0800 1005759
Italy 00800 0225 5734
Spain 00800 0225 5734
United Kingdom 00800 0225 5734
USA & Canada 1 800 228 5166
www.sdi.com.au

EC REP
SDI Germany GmbH
Hansestrasse 85
51149 Cologne
Germany
T+49 2203 9255 0

Last Revision : 02/2021

PAK41398 A